


Product Specifications

| Type | Specifications | Size (Fr) | Effective Length (mm) | Specifications | Size (Fr) | Effective Length (mm) | Specifications | Size (Fr) | Effective Length (mm) | Specifications | Size (Fr) | Effective Length (mm) |
|------|----------------|-----------|-----------------------|----------------|-----------|-----------------------|----------------|-----------|-----------------------|----------------|-----------|-----------------------|
| B | QB040450A | 4 | 450 | QB060900A | 6 | 900 | QB100450A | 10 | 450 | QB060900C | 6 | 900 |
| | QB040550A | 4 | 550 | QB061100A | 6 | 1100 | QB120450A | 12 | 450 | QB070450C | 7 | 450 |
| | QB040700A | 4 | 700 | QB070450A | 7 | 450 | QB060400B | 6 | 400 | QB070550C | 7 | 550 |
| | QB040900A | 4 | 900 | QB070550A | 7 | 550 | QB070400B | 7 | 400 | QB070700C | 7 | 700 |
| | QB041100A | 4 | 1100 | QB070700A | 7 | 700 | QB080400B | 8 | 400 | QB080450C | 8 | 450 |
| | QB050450A | 5 | 450 | QB070800A | 7 | 800 | QB040550C | 4 | 550 | QB080550C | 8 | 550 |
| | QB050550A | 5 | 550 | QB070900A | 7 | 900 | QB040700C | 4 | 700 | QB080700C | 8 | 700 |
| | QB050700A | 5 | 700 | QB071100A | 7 | 1100 | QB050450C | 5 | 450 | QB090450C | 9 | 450 |
| | QB050900A | 5 | 900 | QB080450A | 8 | 450 | QB050550C | 5 | 550 | QB100450C | 10 | 450 |
| | QB051100A | 5 | 1100 | QB080550A | 8 | 550 | QB050700C | 5 | 700 | QB120450C | 12 | 450 |
| | QB060450A | 6 | 450 | QB080700A | 8 | 700 | QB050900C | 5 | 900 | QB050450D | 5 | 450 |
| | QB060550A | 6 | 550 | QB090450A | 9 | 450 | QB060450C | 6 | 450 | QB060450D | 6 | 450 |
| | QB060700A | 6 | 700 | QB090550A | 9 | 550 | QB060550C | 6 | 550 | QB070450D | 7 | 450 |
| | QB060800A | 6 | 800 | QB090700A | 9 | 700 | QB060700C | 6 | 700 | — | — | — |
| C | QC040450A | 4 | 450 | QC060900A | 6 | 900 | QC100450A | 10 | 450 | QC060900C | 6 | 900 |
| | QC040550A | 4 | 550 | QC061100A | 6 | 1100 | QC120450A | 12 | 450 | QC070450C | 7 | 450 |
| | QC040700A | 4 | 700 | QC070450A | 7 | 450 | QC060400B | 6 | 400 | QC070550C | 7 | 550 |
| | QC040900A | 4 | 900 | QC070550A | 7 | 550 | QC070400B | 7 | 400 | QC070700C | 7 | 700 |
| | QC041100A | 4 | 1100 | QC070700A | 7 | 700 | QC080400B | 8 | 400 | QC080450C | 8 | 450 |
| | QC050450A | 5 | 450 | QC070800A | 7 | 800 | QC040550C | 4 | 550 | QC080550C | 8 | 550 |
| | QC050550A | 5 | 550 | QC070900A | 7 | 900 | QC040700C | 4 | 700 | QC080700C | 8 | 700 |
| | QC050700A | 5 | 700 | QC071100A | 7 | 1100 | QC050450C | 5 | 450 | QC090450C | 9 | 450 |
| | QC050900A | 5 | 900 | QC080450A | 8 | 450 | QC050550C | 5 | 550 | QC100450C | 10 | 450 |
| | QC051100A | 5 | 1100 | QC080550A | 8 | 550 | QC050700C | 5 | 700 | QC120450C | 12 | 450 |
| | QC060450A | 6 | 450 | QC080700A | 8 | 700 | QC050900C | 5 | 900 | QC050450D | 5 | 450 |
| | QC060550A | 6 | 550 | QC090450A | 9 | 450 | QC060450C | 6 | 450 | QC060450D | 6 | 450 |
| | QC060700A | 6 | 700 | QC090550A | 9 | 550 | QC060550C | 6 | 550 | QC070450D | 7 | 450 |
| | QC060800A | 6 | 800 | QC090700A | 9 | 700 | QC060700C | 6 | 700 | — | — | — |
| D | QD040900A | 4 | 900 | QD060900A | 6 | 900 | QD080900A | 8 | 900 | QD070550C | 7 | 550 |
| | QD041100A | 4 | 1100 | QD061100A | 6 | 1100 | QD060450B | 6 | 450 | QD050550D | 5 | 550 |
| | QD050900A | 5 | 900 | QD070800A | 7 | 800 | QD070450B | 7 | 450 | QD060550D | 6 | 550 |
| | QD051100A | 5 | 1100 | QD070900A | 7 | 900 | QD050550C | 5 | 550 | QD070550D | 7 | 550 |
| | QD060800A | 6 | 800 | QD071100A | 7 | 1100 | QD060550C | 6 | 550 | — | — | — |

Four tip shapes of Catheter sheath sets Type B/C/D

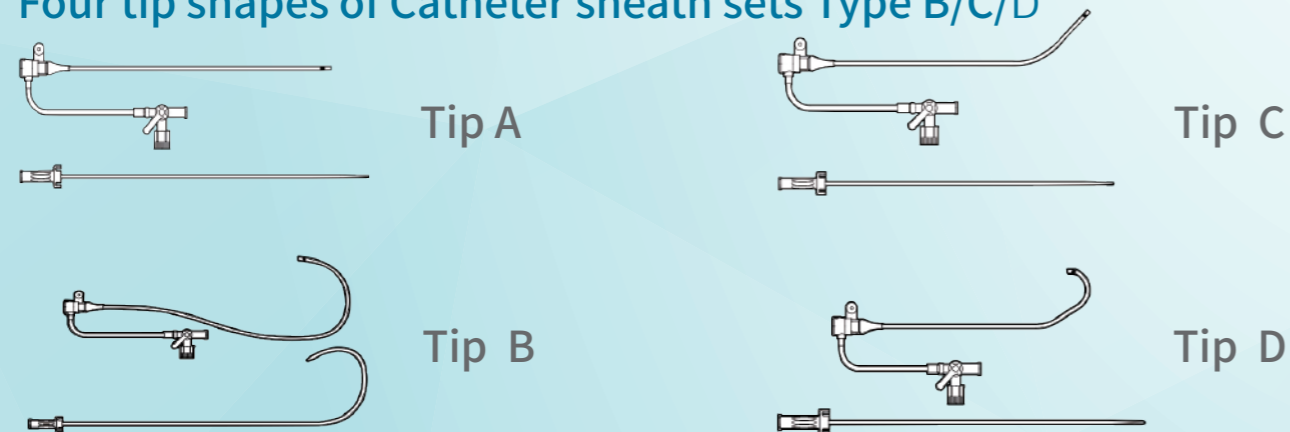


Zhejiang Barty Medical Technology Co.,Ltd.

Add : South 2F, Building 2-1, No.20 Street, ETDZ, Hangzhou, Zhejiang, China
 Tel : 0571-883737191
 Fax : 0571-88373716
 www.vascularbarty.com Web : //www.bartymedical.com
 Email : info@bartymedical.com



HEETY™

Catheter Sheath Sets Type B/C/D

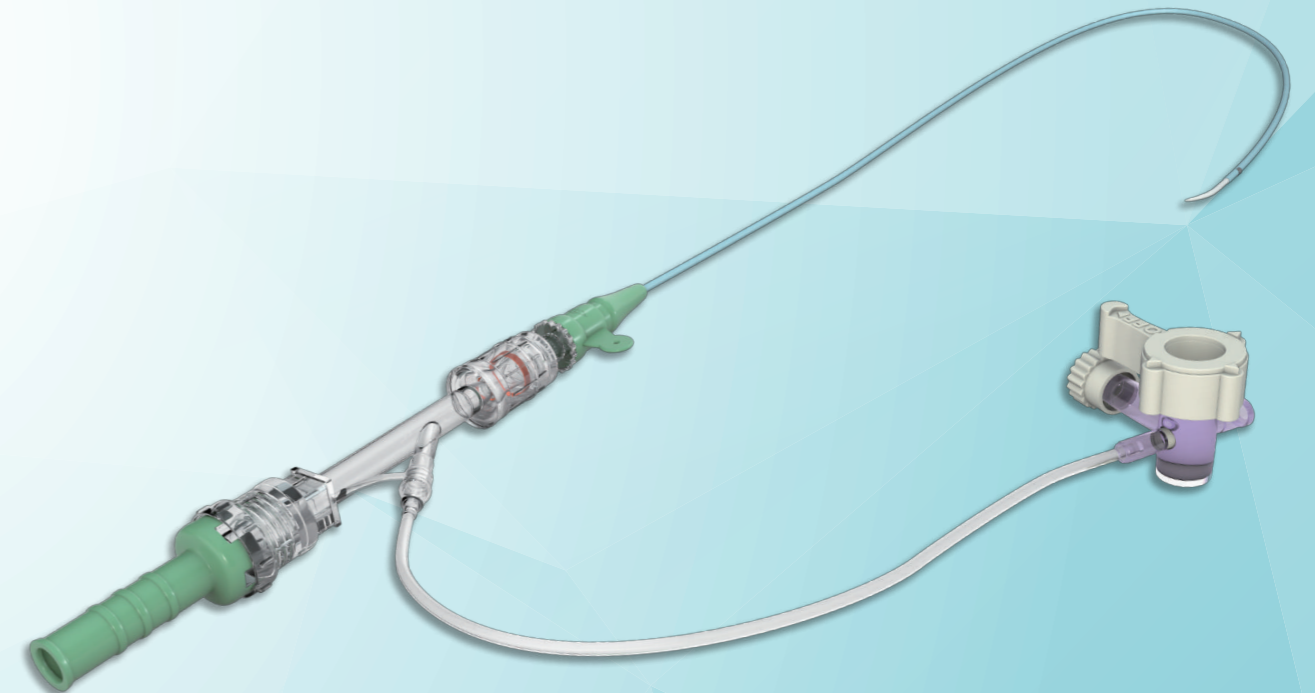

Distributor in Czech Republic:

Endolab vascular s.r.o.

Add: Sokolovská, 428/130, Praha 8 - Karlín, 186 00, Czech Republic

Tel: +420 722 501 700, +420 723 833 454


Web: www.endolaba.com

Email: endolaba@gmail.com

HEETY™

Catheter Sheath Sets Type B/C/D

Variety Of Specifications And Types For More Applications

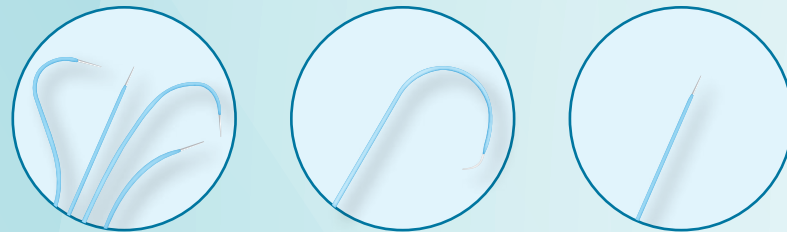


Type D: Y-connector

Design of Y-valve prevents aeroembolism, reduces blood leak, being convenient for operation. No Y-connector in Type A/B/C.

Excellent Flexural-resistance And compliance

Excellent flexural-resistance and compliance:
PTFE inner layer ensures smooth advancement;
Stainless steel inner layer ensures excellent flexural-resistance;
PA outer layer ensures smooth advancement inside blood vessel.

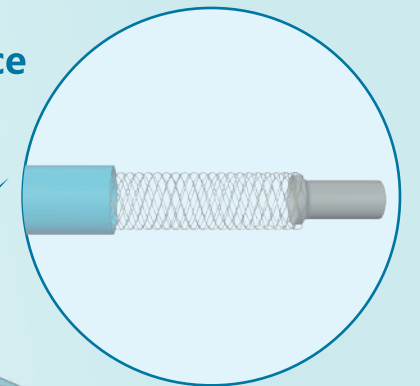


Hydrophilic Coating For Easy Insertion

Hydrophilic coating over sheath surface to provide a smooth surface for insertion.

Safe And Humanized Design To Lower Operational Risk

Perfect matching between sheath and dilator to reduce damage to blood vessel during operation;
Mark band with non transparent X ray on sheath distal end to locate.

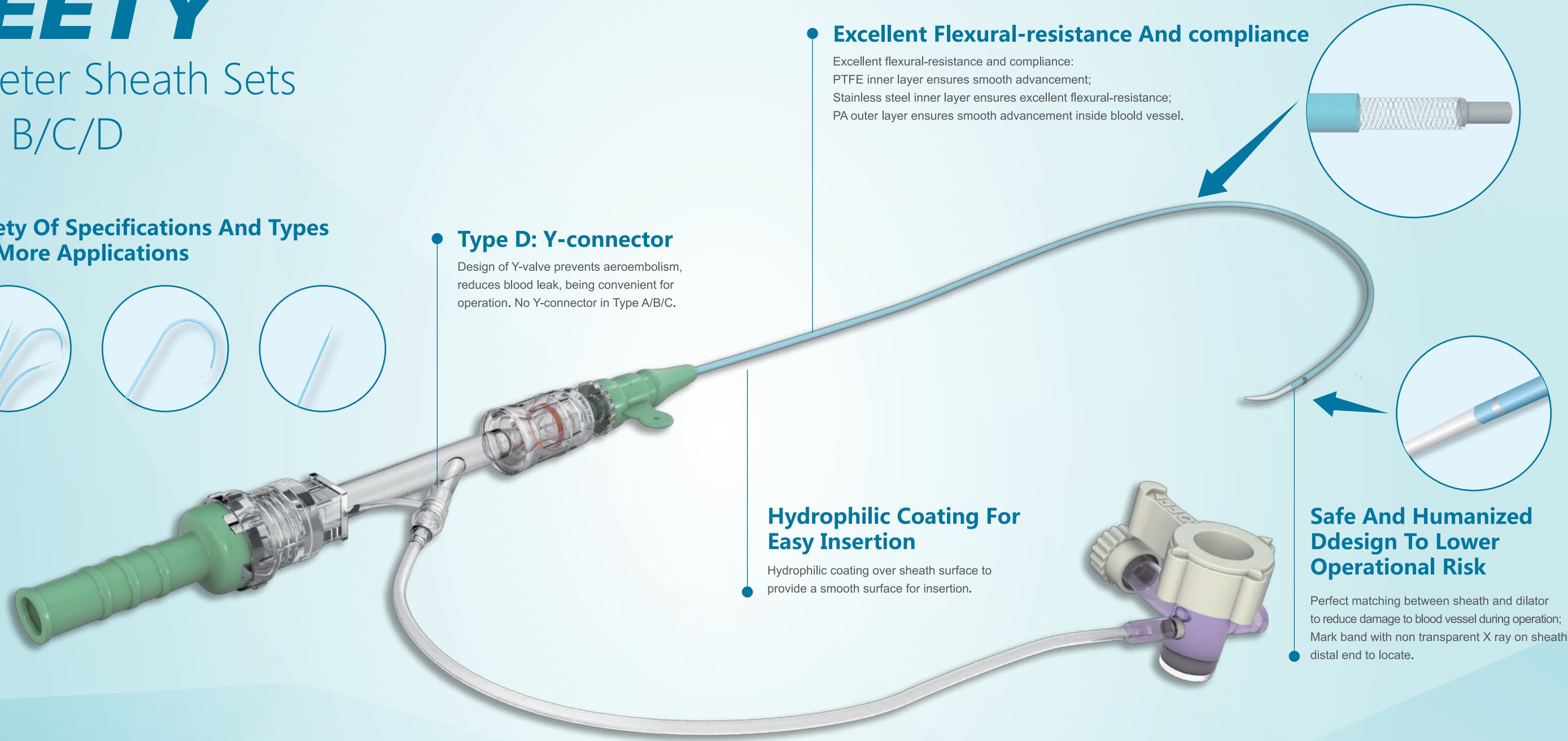

| Type | Product Configuration | | | | | |
|------|-------------------------------|---------|-------------------|---------------|------------|---------|
| B | Catheter Sheath | Dilator | — | — | — | — |
| C | Catheter Sheath | Dilator | Introducer Needle | Scalpel blade | Guide Wire | Syringe |
| D | Catheter Sheath (Y-connector) | Dilator | — | — | — | — |

Application

This product is intended to percutaneously insert into blood vessel during interventional operation, providing channel for angiographic catheter, guide catheter, balloon dilatation catheter and other similar devices into artery and vein.